

FRED® EASYPORT® PLUS


Small and powerful defibrillator for anyone - anytime - anywhere



FRED® EASYPORT® PLUS

The ideal companion on the go

FRED® easyport® plus, the pocket-size defibrillator, fits easily into every doctor's bag and emergency kit. It is the ideal companion during any intervention.

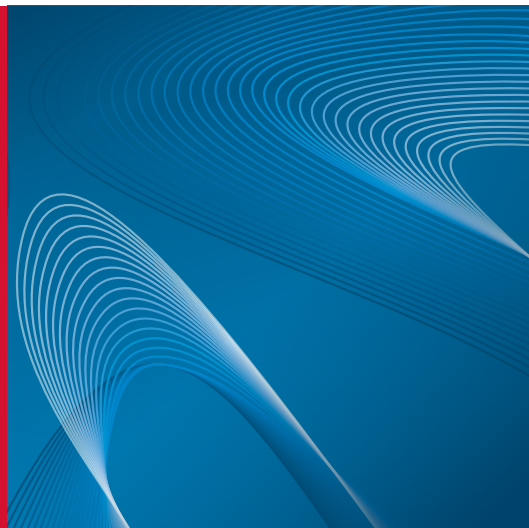

IT CONCERNS ALL OF US

Every year, more than 7 million people fall victim to a sudden cardiac arrest worldwide. A sudden cardiac arrest can hit anyone, anywhere and at any time.

Defibrillation delivers electrical energy to the heart, restoring a proper heart rhythm. A defibrillator in combination with well-performed chest compressions is the only therapy that can save a patient's life.

EASY TO CARRY

As a pocket-size defibrillator, the FRED easyport plus weighs only 780 grams and combines small size with great functionalities.

SHOCK AS SOON AS IT IS REQUIRED

International guidelines emphasise the importance of minimising the delay in shock delivery.

To reduce this delay, the FRED easyport plus has a quick analysis system and loads the required energy systematically so that the shock can be delivered as soon as the device recommends it.

ADAPTED TO ANY USER

The FRED easyport plus is suitable for frequent as well as non frequent users. It can be operated with a rechargeable or a disposable battery. The external charger can charge two batteries simultaneously.

CPR FEEDBACK

Defibrillation alone is not enough. To maintain a good supply of oxygen to the organs, cardiopulmonary resuscitation (CPR) needs to be performed correctly to restore a proper blood flow.

The FRED easyport plus helps you to achieve highest efficiency in chest compressions by:

- indicating the correct depth, frequency and recoil*
- setting a regular rhythm through the metronome

When using the ARGUS LifePoint sensor, the rescuer benefits from visual assistance through on-screen colour-coded messages while performing chest compressions.

*Recoil not available in France, Germany, UK and USA



INFORMATION TRACKING

The FRED easyport plus records intervention data, allowing doctors and emergency physicians to get a better overview of the patient's situation.

- ❖ CPR quality
- ❖ ECG curve
- ❖ Impedance curve
- ❖ Events

INFORMATION RETRIEVAL

Thanks to USB and Bluetooth functionalities, saved files can easily be reviewed and managed by SCHILLER's data management solution SEMA.

FIND THE DEVICE THAT SUITS YOU

The FRED easyport plus is available in three versions to fit the different user levels.



AED FIRST

The FRED easyport plus semi-automatic AED with its grey cover requires the rescuer's confirmation to deliver a shock.

An ECG curve can be displayed on the large LCD colour screen.

Suitable for first responders.



AED AUTO

The fully automatic FRED easyport plus does not require any action from the rescuer. After the heart rate analysis, the device determines whether a shock is necessary. If it is, the shock is delivered automatically.

An ECG curve can be displayed on the large LCD colour screen.

Suitable for trained lay responders.



AED MANUAL

Distinguishable by its red cover, the manual FRED easyport plus allows doctors and paramedics to switch to manual override. Thanks to the ECG display, healthcare professionals can determine if a shock is necessary.

Suitable for physicians.



SCHILLER
FRED easyport plus

00:00:04

PLUG AND APPLY ELECTRODES

1

2

3

ELECTRODES
FRED PA-1,
RFID
ADULT

1

2

3

4

SMALL AND POWERFUL

The FRED easyport plus weighs only 780 grams and delivers shocks of up to 200 joules. The Multipulse Biowave® technology ensures an effective and safe defibrillation shock while minimising damage to the myocardium.

GAIN PRECIOUS SECONDS

The defibrillation electrodes are pre-connected for a quick application on the patient's chest. Heart rate analysis starts right away to reduce time until the shock is delivered.

ALWAYS READY

The FRED easyport plus performs regular automatic self-tests to ensure it works correctly. Its status is indicated by the flashing on/off button.

DATA TRANSMISSION TO SEMA

Intervention data can be sent via USB or Bluetooth for review by healthcare professionals.

DESIGNED TO LAST

With IP44 classification, the FRED easyport plus resists entry of splashing water.





POCKET-SIZE DEFIBRILLATOR

The FRED easyport plus measures only 150 x 143 millimeters – this photo shows its original size (1:1)!

SUITABLE FOR ALL ENVIRONMENTS

The FRED easyport plus has a large colour LCD screen which automatically adjusts its brightness to help visualise indications at a glance.

INTUITIVE PAEDIATRIC MODE

While using the same electrodes, just press this button to change the patient mode.

CPR FEEDBACK

The integrated metronome sets the correct rhythm and guides the rescuer during the intervention.

Additionally, the ARGUS LifePoint sensor provides accurate visual assistance through on-screen colour-coded messages to ensure most effective chest compressions.

OPTIMISED USER INTERFACE

Three simple steps with illustrations and spoken instructions for a clear understanding.



Americas
SCHILLER Americas Inc.
 Doral, Florida 33172

North America:
 Phone +1 786 845 06 20
 Fax +1 786 845 06 02
 sales@schilleramericas.com
 www.schilleramericas.com

Latin America & Caribbean
 Phone +1 305 591 11 21
 Fax +1 786 845 06 02
 sales@schilleramericas.com
 www.schilleramericas.com



Asia
SCHILLER Asia-Pacific / Malaysia
 52200 Kuala Lumpur, Malaysia
 Phone +603 6272 3033
 Fax +603 6272 2030
 sales@schiller.com.my
 www.schiller-asia.com



Austria
SCHILLER Handelsgesellschaft m.b.H.
 A-4040 Linz
 Phone +43 732 70 99 0
 Fax +43 732 757 000
 sales@schiller.at
 www.schiller.at



China
Alfred Schiller (Beijing) Medical Equipment Co.,Ltd.
 100015 Beijing, China
 Phone +86 010 52007020
 Fax +86 010 52007020-8016
 info@schillerchina.com
 www.schillermedical.cn



Croatia
SCHILLER d.o.o.
 10000 Zagreb
 Phone +385 1 309 66 59
 Fax +385 1 309 66 60
 info@schillerzg.hr
 www.schiller.hr



France
SCHILLER Médical
 F-67160 Wissembourg
 Phone +33 3 88 63 36 00
 Fax +33 3 88 63 36 49
 info@schiller.fr
 www.schiller-medical.com



France (distribution France)
SCHILLER France S.A.S.
 F-77608 Bussy St Georges
 Phone +33 1 64 66 50 00
 Fax +33 1 64 66 50 10
 contact@schillerfrance.fr
 www.schiller-france.com



Germany
SCHILLER Medizintechnik GmbH
 D-85622 Feldkirchen b. München
 Phone +49 89 62 99 81 0
 Fax +49 89 609 50 90
 info@schillermed.de
 www.schillermed.de



India
SCHILLER Healthcare India Pvt. Ltd.
 Mumbai - 400 059, India
 Phone +91 22 6152 3333 / 2920 9141
 Fax +91 22 2920 9142
 sales@schillerindia.com
 www.schillerindia.com



Poland
SCHILLER Poland Sp. z o.o.
 PL-02-729 Warszawa
 Phone +48 22 843 20 89 / +48 22 647 35 90
 Fax +48 22 843 20 89
 schiller@schiller.pl
 www.schiller.pl



Russia & C.I.S.
AO SCHILLER.RU
 125040 Moscow, Russia
 Phone +7 (495) 970 11 33
 Fax +7 (495) 956 29 10
 mail@schiller.ru
 www.schiller.ru



Serbia
SCHILLER d.o.o.
 11010 Beograd
 Phone +381 11 39 79 508
 Fax +381 11 39 79 518
 info@schiller.rs
 www.schiller.rs



Slovenia
SCHILLER d.o.o.
 2310 Slovenska Bistrica
 Phone +386 2 843 00 56
 Fax +386 2 843 00 57
 info@schiller.si
 www.schiller.si



Spain
SCHILLER ESPAÑA, S.A.
 ES-28232 Las Rozas/Madrid
 Phone +34 91 713 01 76
 Fax +34 91 355 79 33
 schiller@schiller.es
 www.schiller.es



Switzerland
SCHILLER-Reomed AG
 CH-8912 Obfelden
 Phone +41 44 744 30 00
 Fax +41 44 740 37 10
 sales@schiller-reomed.ch
 www.schiller-reomed.ch



Turkey
SCHILLER TÜRKIYE
 Okmeydani-Sisli - Istanbul
 Phone +90 212 210 8681 (pbx)
 Fax +90 212 210 8684
 info@schiller.com.tr
 www.schiller-turkiye.com



United Kingdom
SCHILLER UK
 Bellshill, ML4 3PR
 Phone +44 1698 744 505
 Fax +44 1698 744 474
 sales@schilleruk.com
 www.schilleruk.com

Device availability in your market is subject to regulatory approval.

